



 **HYGENIC[®]**
DENTAL DAM

**Many Paths Lead to Success –
Take the Dry One!**

Improved clinical results with greater moisture and infection control



Excellent visibility and protection during root canal preparation



Patient optimally protected against aggressive rinsing agents



Use rubber dam even for subgingival preparations

Save time

The few minutes that it takes to apply dental dam is more than compensated for by the time saved during treatment. Properly placed dam improves visibility and protects cheeks, lips, and tongue. Improved access and visibility allows treatment to be completed in a timely and efficient manner.

Optimum protection for patient...

Your patients are protected against aspiration or ingestion of foreign objects. Especially in endodontic treatment or esthetic restorations, it is essential to establish a dry and clean operating field - a basis for durable clinical success.


...and dental staff

Protect yourself and your dental team against viruses and bacteria which become airborne when using rotary instruments. Dental dam limits the risk of infection from airborne saliva and blood particles.

Quick and easy placement...



Stretch dental dam over a dam frame and punch out the pre-marked holes.

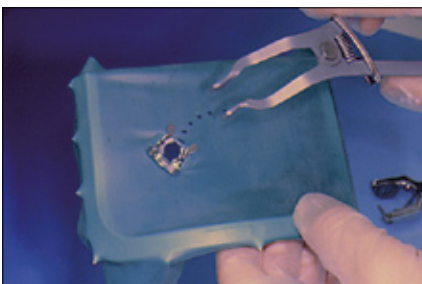

Insert winged clamp and place with clamp forceps.



Position clamp with dental dam frame in one easy step.



An alternative to clamps: Wedjets® stabilizing cord provide secure and gentle retention.

HYGENIC® – A diverse and complete line of Dental Dam and accessories



HYGENIC® Dental Dam*

High quality, pure natural rubber latex - strong and tear resistant.

HYGENIC® offers a comprehensive dental dam line:

Several sizes of pre-cut dam, colors and thicknesses are available to suit your individual needs. Rubber dam in rolls can be cut to any individual size.

Fiesta® Dental Dam*

with fruit scent for even more patient comfort is available in assorted pastel colors.

Double security for persons suffering from allergies

HYGENIC® Non-Latex Dental Dam

Latex allergies have become increasingly common among patients and dental staff. Special solutions are required. Two popular non-latex dental dam choices are available:

HYGENIC® Non-Latex Dental Dam is offered in both 127 × 127 mm and 152 × 152 mm.

Non-Latex Flexi Dam® offers greater strength. Available in 152 × 152 mm Purple or Tealgreen.

Framed Non-Latex Flexi Dam®

The Hygenic Non-Latex Flexi Dam is available with an ultra-convenient, built-in-frame. The flexible frame is designed with a convenient working size of 100 mm x 105 mm to ensure for easy placement without getting in the way. The smooth surface of the plastic frame helps to maximize patient comfort when positioned on their skin.

Latex allergies and sensitizing can be avoided.

Both latex-free alternatives are convincing:

- More tear resistant and elastic than latex dam
- No latex allergy
- Powder-free
- No offensive odor
- Excellent recovery properties

GET ORGANIZED... with Fiesta® Color Coded Clamps

HYGENIC® Fiesta®

Color Coded Clamps

Each clamp is individually color coded for unmistakable identity. Experienced and novice assistants will appreciate this time saving identification when using Fiesta Color Coded Clamps. For convenient storage, the 9 color-coded clamps are placed on 9 matching color coded pegs on the stainless steel organizer.



*These products contain natural rubber latex and may cause allergies.



Dental Dam Clamps: Secure with gentle pressure



Winged clamps

The Dam is hinged to the wings of the clamp outside the patient's mouth. The dental dam, dental dam clamp and frame are all simultaneously positioned into the patient's mouth. Once applied, the dam is stretched off the wings of the clamp.



Wingless clamps

First the clamp is positioned intraorally, then the dental dam stretched over the clamp and tooth. Wingless clamps are designated with a „W“ (wingless) for easy identification.

FIESTA Color Coded Matte Finish Winged and Wingless Clamps

ANTERIOR



#9*

#W9**

MOLAR



#W3**

#7*

small lower molars

lower molars

PREMOLAR



#2*

#W2**

large bicuspid

#2A**

#W2A**

bicuspid



#W8**
upper molars

#4*
small upper molars

#W56**
upper & lower molars

SERRATED JAWS

- Serrations for improved stabilization



#12A*

#13A*

lower right molars /
upper left molars

lower left molars /
upper right molars

For irregularly shaped, structurally
compromised or partially erupted molars



#8A*

#W8A**

#14A*

#W14A**

small

large

* Included in **Fiesta 9** Winged Clamp-Pak

** Included in the **Fiesta 9W** Wingless Clamp-Pak

Hygienic Gloss Finish Winged and Wingless Clamps

ANTERIOR



#9*

#W9

#212

MOLAR



#3

#W3

#7*

#W7

small lower molars

lower molars

PREMOLAR



#00

#W00

#1*

upper & lower bicuspid

upper bicuspid



#8

#W8

#4*

#W56

lower molars

small upper molars

upper & lower molars



#2*

#W2

#2A

#W2A

large bicuspid

bicuspid

For irregularly shaped, structurally
compromised or partially erupted molars



#8A*

#W8A

#14

#W14

small

medium

SERRATED JAWS

- Serrations for improved stabilization



#12A*

#13A*

lower right molars /
upper left molars

lower right molars /
upper left molars



#14A*

#W14A

large

* included in **System 9** Clamp-Pak

BRINKER TISSUE RETRACTORS

- for gingival retraction



#B1

Lower Molars

#B2

upper left molars

#B3

upper right molars

#B4

anterior & canines

#B5

class V restorations on all teeth

#B6

HYGENIC® DENTAL DAM

You Value Security – Your Patients Value You!

«Although rubber dam placement may consume some time - the advantages of being able to work in a clean, dry operation field eventually saves me time.

And it makes treatment much more comfortable for my patients.»

Dr. Spalek, Stuttgart



Dental Dam Accessories

Wedjets® – An alternative to conventional clamps

Wedjets® elastic stabilizing cord provide easy and quick retention of the dam. No clamps are required. Available in latex and non-latex.




ORA-SHIELD® – Dental Dam Napkins Additional comfort for patients!

Application of napkins under the dental dam both prevent direct contact of the Dental Dam with the sensitive skin of the patient's face and absorbs moisture.



Dental Dam Stamps or Templates

Accurately mark tooth positions.

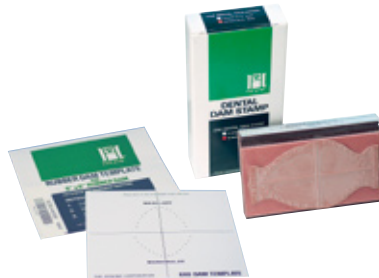

Dam Frames -

Three types to choose from!

- Ostby plastic frame
- U-shaped plastic frame
- U-shaped metal frame

All frames may be autoclaved.



Dental Dam Forceps

Safely position clamps with our light - weight, matte finished forceps



Dental Dam Punch

Ainsworth-style punch with five hole sizes, punches cleanly through any gauge dam



Extended Reach Punch

Hygenic Extended Reach Punch has 12 mm greater reach over the traditional Ainsworth style punch. It is ideal for use with any pre-framed dental dam. The dagger and dies are easy to change without any additional tools. The punch comes with 2 different sized dies: #2 for anterior teeth and #4 for posterior teeth.



*These products contain natural rubber latex and may cause allergies.

HYGENIC®

Ordering information

DENTAL DAM

Dental Dam* ready cut, for children
127 x 127 mm, 52 sheets¹⁾

Thin: H00523 Light Medium: H00524 Light
H02141 Green H00529 Dark
H02142 Green
H03529 Blue

Heavy: H00525 Light Extra heavy: H00526 Light
H00530 Dark H00531 Dark
H02143 Green H02144 Green
H07314 Blue

Dental Dam* ready cut, for adults
152 x 152 mm, 36 sheets¹⁾

Thin: H00533 Light Medium: H00534 Light
H00538 Dark H00539 Dark
H02146 Green H02147 Green
H03530 Blue

Heavy: H00535 Light Extra heavy: H00536 Light
H00540 Dark H00541 Dark
H02148 Green H02149 Green
H07315 Blue

¹⁾Bulkpackages available on request

Roll Dam* for children
127 mm x approx. 6.7 m

Thin: H00503 Light Medium: H00504 Light

Roll Dam* for adults

152 mm x approx. 5.5 m
Light: H02174 Green
Medium: H00514 Light
H00519 Dark
H02175 Green
Extra heavy: H02177 Green

Non-Latex Flexi Dam*

152 x 152 mm, 30 sheets, for adults
Medium: H09945 Purple
H09946 Tealgreen

Framed Non-Latex Flexi Dam*

H00750 purple 20 pc

Non-Latex Dental Dam

152 x 152 mm, 15 sheets, for adults
Medium: H09105 Tealgreen
127 x 127 mm, 15 sheets, for children
Medium: H09928 Tealgreen

Fiesta Dental Dam*, for children

127 x 127 mm, 52 sheets
Thin: H04640 assorted colours
Medium: H04641 assorted colours
Heavy: H04643 assorted colours

Fiesta Dental Dam*, for adults

152 x 152 mm, 36 sheets
Thin: H04639 assorted colours
Medium: H04642 assorted colours
Heavy: H04644 assorted colours

DENTAL DAM KITS

Wingless Complete Kit

H02790 36 sheets HYGENIC® Dental Dam* (green, medium,
152 x 152 mm), 152 mm Plastic Frame (H01414),
Fiesta 9W Wingless Color Coded Clamp Pak,
Template, Wedjets® Cord (Small),
Dental Dam Punch, Dental Dam Forceps

Winged Complete Kit

H02778 36 sheets HYGENIC® Dental Dam* (green, medium,
152 x 152 mm), 152 mm Plastic Frame (H01414),
Fiesta 9 Winged Color Coded Clamp Pak,
Template, Wedjets® Cord (Small),
Dental Dam Punch, Dental Dam Forceps

Winged Complete Kit (packed in Plastic Carry Case)

H02782 36 sheets HYGENIC® Dental Dam* (green, medium,
152 x 152 mm), 152 mm Plastic Frame (H01414),
Fiesta 9 Winged Color Coded Clamp Pak,
Template, Wedjets® Cord (Small), Dental Dam Punch,
Dental Dam Forceps, Plastic Carry Case

Dental Dam Kit with Brinker Clamps (Gloss Finish Clamps)

H09715 52 sheets HYGENIC® Dental Dam* (green, medium,
127 x 127 mm), 127 mm Plastic Frame (H01416),
Fiesta 9W Wingless Color Coded Clamp Pak,
Template, Wedjets®, Punch, Forceps

Wingless Standard Kit (w.o. punch & forceps)

H02792 36 sheets HYGENIC® Dental Dam* (green, medium,
152 x 152 mm), 152 mm Plastic Frame (H01414),
Fiesta 9W Wingless Color Coded Clamp Pak,
Template, Wedjets® Cord (Small)

Wingless Standard Kit (w.o. punch & forceps) – Metal frame

H02793 36 sheets HYGENIC® Dental Dam* (green, medium,
152 x 152 mm), 152 mm Metal Frame (H07083),
Fiesta 9W Wingless Color Coded Clamp Pak,
Template, Wedjets® Cord (Small)

Winged Standard Kit (w.o. punch & forceps)

H02779 36 sheets HYGENIC® Dental Dam* (green, medium,
152 x 152 mm), 152 mm Plastic Frame (H01414),
Fiesta 9 Winged Color Coded Clamp Pak,
Template, Wedjets® Cord (Small)

DENTAL DAM ACCESSORIES

Wedjets 2,1 m

Latex*
H08040 Extra Small, Blue
H06522 Small, Yellow
H06523 Large, Orange

Non-Latex

H06524 Small, Lime Green
H06525 Large, Aqua Blue

Dental Dam Frame

H01416 Plastic, for children (127 mm)
H01414 Plastic, for adults (152 mm)
H00560 Ostby frame, plastic, for adults
H07082 Metal, for children (127 mm)
H07083 Metal, for adults (152 mm)

Dental Dam Forceps

H01252 stainless steel

Dental Dam Punch

H01251 stainless steel

Ora-Shield Dental Dam Napkins, 50 pieces

H01415 fits Frame size 127 mm and 152 mm
H00841 Fits Holder (large size)

Extended Reach Punch

H01253 Punch with dagger and one #2 and #4 die of each
H01254 #2 die 1 pc.
H01255 #4 die 1 pc.
H01256 dagger 1 pc.

Dental Dam Stamp

H01442 for children
H01443 for adults

Dental Dam Template

H01440 for children
H01441 for adults

Clamps Organizer

H09956 Fiesta Color Coded Clamp Organizer with 9 Color
Coded Storage Pegs

DENTAL DAM CLAMPS

(stainless steel)

Clamp-Kits

H02701 System-9-clamp-kit (winged)
Size: 1, 2, 4, 7, 8A, 9, 12A, 13A, 14A

H01054 Brinker Clamp-Set
Size: B1, B2, B3, B4, B5, B6

H09966 Fiesta 9 Winged Color Coded Clamp Pak (Color
Coded Clamp Organizer and Winged sizes: 2, 2A,
4, 7, 8A, 9, 12A, 13A, 14A)

H02801 Fiesta 9W Wingless Color Coded Clamp Pak (Color
Coded Clamp Organizer and Wingless sizes: W2,
W2A, W3, W7, W8, W8A, W9, W14A, W56)

HYGENIC® Dental Dam Clamps, winged

H02762 #00	H02760 #8
H02703 #1	H02759 #8A
H02761 #2	H02763 #9
H05688 #2A	H02706 #12A
H05689 #3	H02705 #13A
H02704 #4	H05690 #14
H02758 #7	H02757 #14A

HYGENIC® Dental Dam Clamps, wingless

H07792 #212	H02767 #W8
H02769 #W00	H02766 #W8A
H02768 #W2	H02770 #W9
H02775 #W2A	H02774 #W14
H03814 #W3	H02764 #W14A
H02765 #W7	H03815 #W56

Brinker Tissue-Retractors-Clamps

H01048 B1	H01051 B4
H01049 B2	H01052 B5
H01050 B3	H01053 B6

Fiesta Color Coded Matte Finish Winged Clamps

H09957 #2	H09960 #7	H09963 #12A
H09958 #2A	H09961 #8A	H09964 #13A
H09959 #4	H09962 #9	H09965 #14A

Fiesta Color Coded Matte Finish Wingless Clamps

H02794 #W2	H02797 #W7	H02800 #W9
H02795 #W2A	H02798 #W8	H02802 #W14A
H02796 #W3	H02799 #W8A	H02803 #W56

*These products contain natural rubber latex and may cause allergies.

www.coltenewhaledent.com
© 2007 Coltène/Whaledent AG

Coltène/Whaledent Inc.

235 Ascot Parkway
Cuyahoga Falls, Ohio 44223/USA
Tel. +1 330 916 8800
Fax +1 330 916 7077
productinfo@coltenewhaledent.com

Coltène/Whaledent Ltd.

The President Suite-A
Kendal House, Burgess Hill
Victoria Way
West Sussex, RH15 9NF/U.K.
Tel. +44 (0)1444 235486
Fax +44 (0)1444 870640
hq@coltenewhaledent.co.uk

Coltène/Whaledent AG

Feldwiesenstrasse 20
9450 Altstätten/Switzerland
Tel. +41 (0)71 757 53 00
Fax +41 (0)71 757 53 01
office@coltenewhaledent.ch

coltène whaledent