

From front...



to back...



*to everywhere
in between.*



The Sof-Lex™ Finishing and Polishing System

3M ESPE



3M ESPE is a recognized leader in dental innovation. Introduced in the mid-1970s, the Sof-Lex™ Finishing and Polishing System was a revolutionary breakthrough – providing dentists with the ultimate advanced restoration finishing solution.

With a color-coded, sequential system of discs, strips and easy-to-use mandrels, the Sof-Lex system continues to provide unparalleled, efficiently achieved results.



Now, the renaissance continues...
by pushing innovation to
the tip with the Sof-Lex
Finishing Brush.

Introducing the Sof-Lex™ Finishing Brush

The Sof-Lex Finishing Brush combines the once separate finishing and polishing procedures of posterior composite restorations into *one, easy step*. No stopping to switch shapes or grits saves time and effort – making restorative procedures much more efficient for you.

The brush bristles contain aluminum oxide abrasives for highly effective polishing. The soft, elastomer bristles are gentle to gingival tissue and conform as they travel over concave and convex surface anatomies to deliver a uniform finish comparable to Sof-Lex discs.

The Sof-Lex Finishing Brush is also *autoclavable*. As part of the 3M ESPE Sof-Lex Finishing and Polishing System, the finishing brush delivers the superior results you demand – in *less time* and with *less effort* than you expect.

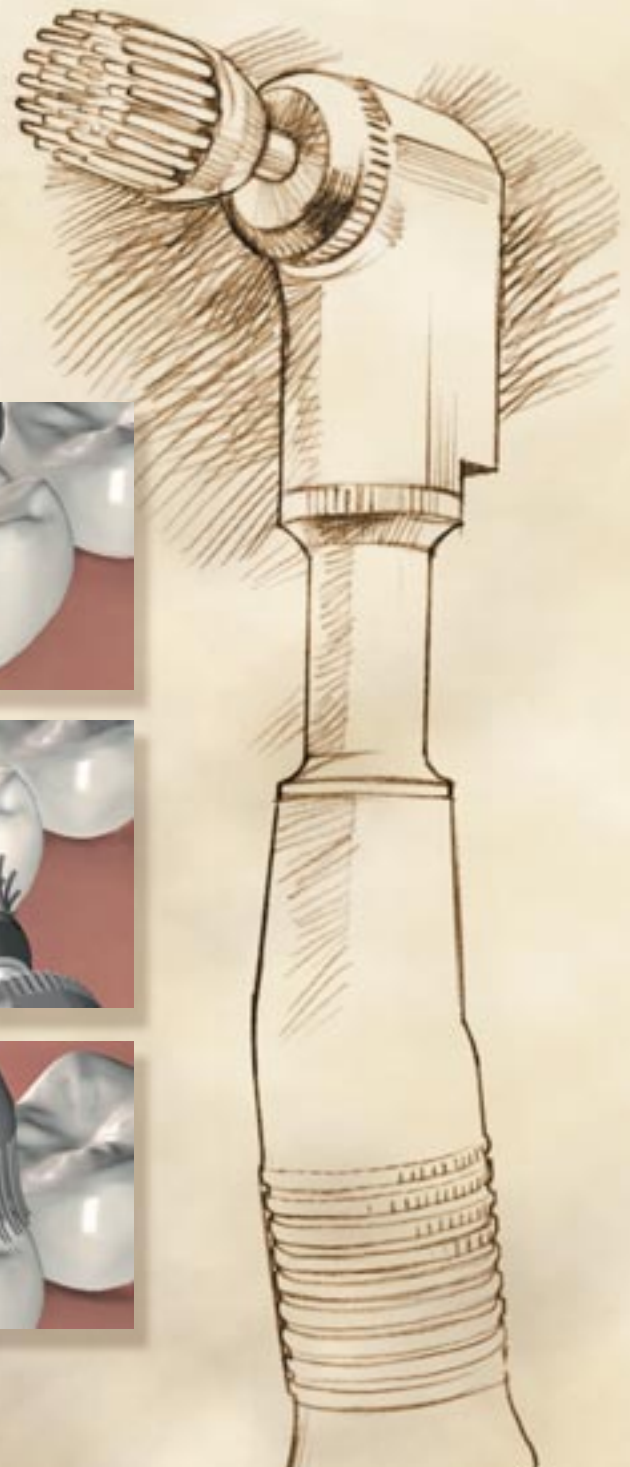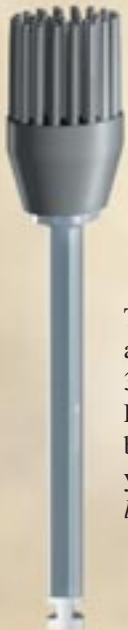

Sof-Lex™ Finishing and Polishing Discs and Strips

Pick up the pace without sacrificing great results.

With coarse to superfine grits clearly color-coded for ease of procedure sequencing and quick, pop-on/off disc interchangeability, Sof-Lex finishing and polishing strips and discs deliver exceptional, *time-proven* results with *time-saving* advantages.



Sof-Lex Discs

- **Reversible Discs**
Makes aligning disc easier for polishing various surfaces.
- **Fast and Easy Disc Change**
Simply snaps securely on or off mandrel with no need for alignment.
- **Color-Coded**
Sequenced from dark (coarse) to light (superfine) for easily followed step-by-step process.
- **Two Sizes**
For a variety of needs.
13mm (1/2 inch) or 9mm (3/8 inch)
- **Maximized Abrasive Surface**
Mandrel is low profile. Disc has a small round eyelet.
- **Choice of disc Thickness and Flexibility**
Applicable to a variety of uses with regular or extra thin construction.



Sof-Lex™ Contouring and Polishing Discs



Sof-Lex™ Extra Thin Contouring and Polishing Discs



Sof-Lex Strips

For finishing and polishing interproximal areas, Sof-Lex strips round-out the system.

- **Easy to Insert**
Abrasive is center gapped on the strip for easy interproximal insertion *without* abrading contacts.
- **Two Grits per Strip**
Each strip performs two functions to save time.
- **Color-Coded**
Sequenced from dark to light for easily followed step-by-step process.
- **Choice of Grits and Widths**
A choice of four grits and two widths.
- **Resists Tearing**
Polyester backing resists tearing while retaining flexibility and being gentle on gingiva.

Sof-Lex™ Finishing and Polishing Strips



Restorative and Adhesive Products

Prompt™ L-Pop™ All-In-One Adhesive

- New formulation – allows the popular adhesive to be cured with any light source – halogen, laser, plasma or LED.
- Yellow tint allows better visibility – offers increased visibility and control while wetting the brush and applying the adhesive to the tooth. Also acts as a visual cue that the product has been activated.
- New extra-fine brush – extra-small (XS) brush with yellow handle available for very small preparations.



Single Bond Dental Adhesive System

- Excellent bonds to both dentin and enamel – thus reducing microleakage.
- Fast, easy technique – saves application time and steps.
- Unique bottle design – easy to open, easy to dispense, and you can see how much is left inside.
- Appropriate for direct and indirect restorations – compatible with 3M ESPE RelyX™ ARC adhesive resin cement making indirect restorations easy, too.



Filtek™ Z250 Universal Restorative

- Excellent esthetics – available in 15 shades to meet your needs.
- Speed and convenience – easy to use. Cures in just 20 seconds* and finishes to a smooth surface.
- Outstanding versatility – useful for any area of the mouth.
- Low polymerization shrinkage – less stress on the tooth, less sensitivity potential for the patient.



Restorative and Adhesive Products

Filtek™ P60 Posterior Restorative

- Ideal handling – excellent packability and resists slumping. Suitable for all Class I and II restorations.
- Outstanding physical properties – excellent wear resistance and overall strength.
- Fast placement technique – takes only 20 seconds to cure an increment of 2.5mm.
- Does not stick to dental instruments – saves time struggling with pull backs.



Filtek™ Flow Flowable Restorative

- “Smart” handling – flows under pressure, but holds its shape and stays in place prior to light cure.
- Ideal flow – viscosity suitable for small Class I, III, IV and shallow Class V restorations.
- Excellent strength – great compressive and tensile strength. Optimal wear resistance.
- Easy to use – just a 20-second cure time and direct syringe delivery.



Elipar™ FreeLight Curing Light

- Cord free – take it from room to room. The base provides optimum recharging, protecting the battery and preventing memory effect.
- Bulb free – no halogen bulbs. No cooling fan. No noise.
- Bulk free – lightweight, ergonomic, easy to handle.
- Effective polymerization – blue LEDs generate a light that closely matches the photoinitiator systems most frequently used in dental materials.



Sof-Lex™ Contouring and Polishing Discs



Sof-Lex™ Contouring and Polishing Discs Kit #1980

- 30 ea. Coarse, Medium, Fine and Superfine Discs – 3/8"
- 30 ea. Coarse, Medium, Fine and Superfine Discs – 1/2"
- 1 RA Mandrel


Sof-Lex™ Extra Thin (XT) Contouring and Polishing Discs



Sof-Lex™ Extra Thin Contouring and Polishing Discs Kit #2380

- 30 ea. Extra Thin Coarse, Medium, Fine and Superfine Discs – 3/8"
- 30 ea. Extra Thin Coarse, Medium, Fine and Superfine Discs – 1/2"
- 1 RA Mandrel

Sof-Lex™ Finishing and Polishing Strips

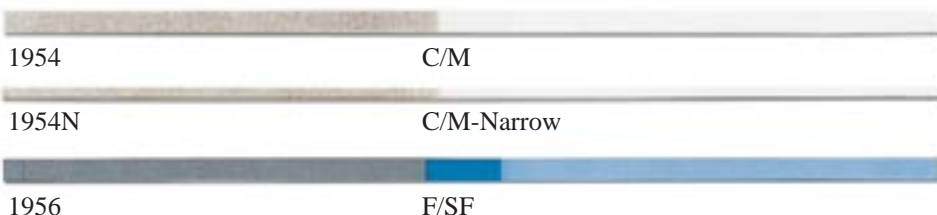

Sof-Lex™ Finishing and Polishing System Kit #2385P

- 120 Sof-Lex Contouring and Polishing Discs
- 120 Sof-Lex Extra Thin Contouring and Polishing Discs
- 45 Sof-Lex Finishing and Polishing Strips
- 1 RA Mandrel

Sof-Lex™ Finishing Brush

Intro Kit #1750

- 12 Sof-Lex Finishing Brushes
- 2 Sof-Lex Finishing Brush Mandrels – Right Angle
- 30 ea. Assorted Sof-Lex Extra Thin (XT) Contouring and Polishing Discs – Coarse, Medium, Fine and Superfine – 1/2"
- 1 Sof-Lex Disc Mandrel – Right Angle



Brush Refill #1751

- 24 Sof-Lex Finishing Brushes



Mandrel Refill #1755 – Right Angle #1756 – Friction Grip

- 3 Sof-Lex Finishing Brush Mandrels

3M ESPE Technical Hotline: 1-800-634-2249

3M ESPE Website: www.3MESPE.com

*Filtek Z250 Universal Restorative shades UD, C4 and B0.5 require a 30-second cure time.

3M ESPE

Dental Products

3M Center
Building 275-2SE-03
St. Paul, MN USA 55144-1000

Dental Products

3M Canada, Inc.

P.O. Box 5757
London, Ontario, Canada N6A 4T1
1-800-265-1840 ext. 6229



10% Post-consumer waste paper

©3M 2002. 3M, ESPE, Sof-Lex, Elipar, Filtek, Scotchbond, Prompt and L-Pop are trademarks of 3M or 3M ESPE AG. 20-2009-3409-2